ageLOC[®]
WellSpa iO[™]

USER'S MANUAL AND DECLARATION OF CONFORMITY

CONTENTS

EN	USER'S MANUAL AND DECLARATION OF CONFORMITY	4
CS	UŽIVATELSKÁ PŘÍRUČKA A PROHLÁŠENÍ O SHODĚ	20
DA	BRUGERMANUAL OG OVERENSSTEMMELSESERKLÆRING	36
DE	GEBRAUCHSANLEITUNG UND KONFORMITÄTSEKTLÄRUNG	52
ES	MANUAL DEL USUARIO Y DECLARACIÓN DE CONFORMIDAD	68
FI	KÄYTTÖOHJEET JA VAATIMUSTENMUKAISUUSVAKUUTUS	84
FR	MANUEL D'UTILISATION ET DÉCLARATION DE CONFORMITÉ	100
HU	HASZNÁLATI ÚTMUTATÓ ÉS MEGFELELŐSÉGI NYILATKOZAT	116
IT	MANUALE UTENTE E DICHIARAZIONE DI CONFORMITÀ	132
NL	GEBRUIKERSHANDLEIDING EN CONFORMITEITSVERKLARING	148
NO	BRUKERVEILEDNING OG SAMSVARERKLÆRING	164
PL	PODRĘCZNIK UŻYTKOWNIKA I DEKLARACJA ZGODNOŚCI	180
PT	MANUAL DO UTILIZADOR E DECLARAÇÃO DE CONFORMIDADE	196
RO	MANUAL DE UTILIZARE ŞI DECLARAŢIE DE CONFORMITATE	212
SK	NÁVOD NA POUŽITIE A VYHLÁSENIE O ZHODE	228
SV	ANVÄNDARMANUAL OCH DEKLARATION OM ÖVERENSSTÄMMELSE	244



TABLE OF CONTENTS

System Contents and Components	5–6	Maintenance	17
Safety Instructions, Warnings, and Precautions	7	Replacement and Warranty Information	17
Charging Your ageLOC® WellSpa iO™	10	Patents	17
Connecting Your ageLOC WellSpa iO to the Nu Skin Vera® App	11	Technical and Regulatory Information	18
Pairing Your ageLOC WellSpa iO	11	Definition of Symbols	18
Using Your ageLOC WellSpa iO	12	Declarations of Conformity	19
Cleaning and Caring for Your ageLOC WellSpa iO	15		
Troubleshooting	15		

SYSTEM CONTENTS


ageLOC® WellSpa iO™ device

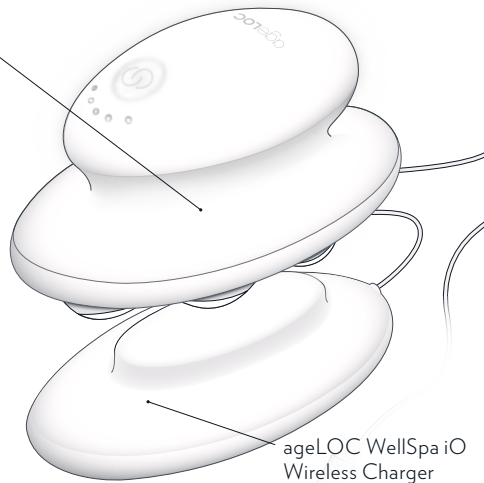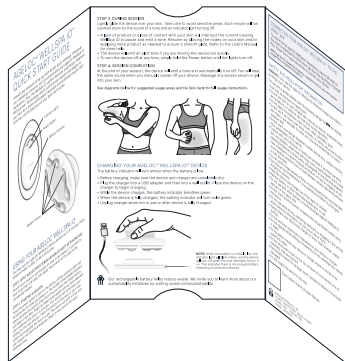Quick Start Guide

USB charging adapter not included.

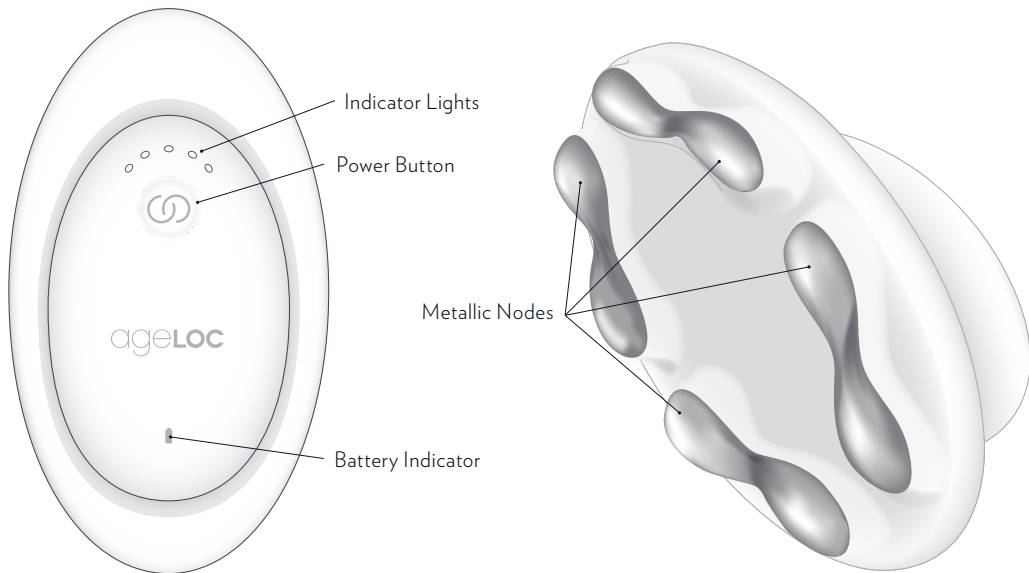SYSTEM COMPONENTS

SAFETY INSTRUCTIONS AND WARNINGS

- This device is for use by adults above 18 years old. It is not intended for use by persons with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. Cleaning and user maintenance shall not be made by children.
- ageLOC® WellSpa iO™ and the ageLOC WellSpa iO Wireless Charger are tested at pressures matching 4000m in altitude.
- ageLOC WellSpa iO contains a lithium-ion battery that is not replaceable. To prevent the possibility of injury, never expose the device to heat. Do not store near a heat source, such as a radiator, fire, or heat vent. Do not leave in a hot vehicle or in an extremely cold environment for prolonged periods of time. Check with your carrier before shipping or flying with this device.
- **WARNING:**  The ageLOC WellSpa iO device is charged using inductive charging. Only use the ageLOC WellSpa iO Wireless Charger provided to recharge your device and a good-quality USB-enable power adapter with the following ratings: 5V  min 500 mA. No other charging accessories should be used with the device.
- Always place your ageLOC WellSpa iO device and ageLOC WellSpa iO Wireless Charger on a dry, heat-resistant, stable, flat surface, especially when charging.
- Do not keep the charger plugged in when not in use or after device has reached full charge.
- Ideal device storage and charging is to occur at room temperature, between 20°C and 24°C (68°F to 75°F).
- Only store the device in environments between 5°C and 42°C (41°F to 107°F).
- Only charge the device in environments between 10°C to 35°C (50°F to 95°F).

SAFETY INSTRUCTIONS AND WARNINGS (Cont.)

- Storage and charging outside of these ranges will adversely affect the battery, metallic nodes, and operation of the device.
- Very high temperatures or hot environments, (e.g., above 50°C/122°F, direct sunlight, within vehicles in extremely hot conditions, etc.) can result in overheating and other catastrophic conditions that severely impact the performance and life of the device.
- Very low temperatures or freezing environments can result in permanent damage to the metallic nodes and operation.

PRECAUTIONS

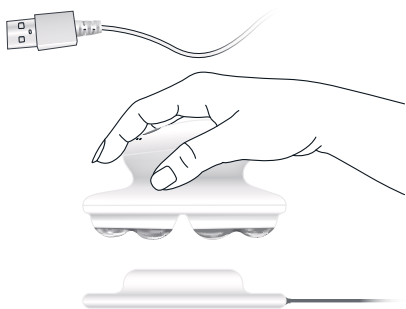
- Keep the device and charger out of the reach of children.
- ageLOC WellSpa iO is for individual use only.
- Use only as directed on healthy skin.
- Periodically inspect for damage; never operate the device if it appears to be damaged or broken to avoid injuries. In the case of any device defect or if the device gets damaged, contact Nu Skin Support Services.
- Do not use over or near irritated skin, cancerous lesions, diseased tissue, skin eruptions, or areas of skin that lack normal sensation as it could result in irritation.
- Do not apply to head/face, neck, throat, chest, décolleté, groin.
- Do not use if you have open sores, wounds, or rashes, are prone to pityriasis rosea, have allergies to metal, or if you have problematic or overly sensitive skin.

- Do not use over swollen, red, infected, or inflamed areas; or over skin eruptions (e.g., phlebitis, thrombophlebitis, varicose veins).
- Do not use if you have an implanted cardiac device (e.g., cardiac monitors, ICDs, pacemakers, ECG alarms).
- Consult with your physician or healthcare provider before using the device if you have sensitive skin, are epileptic, wear body piercings, have a metal implant, or are ill, have bleeding tendencies, bruise easily, are pregnant or may become pregnant, have undergone a recent surgical procedure, or have an ongoing or previous medical condition such as heart disease, stroke, high blood pressure or liver disorders.
- Do not use for longer periods of time or more frequently than recommended.
- Discontinue use immediately if you experience irritation or prolonged stinging or redness and seek medical advice if necessary.
- Remove jewellery or other metal adornments from the surface of the skin prior to use.
- Do not use with other devices, including multiple ageLOC WellSpa iO devices, at the same time.
- Do not touch another person or electrically live parts/items while using the device.
- Do not use while driving, operating machinery, or engaging in any activity where electrical stimulation can put you at risk of injury.
- Do not use harsh chemicals or abrasives on your ageLOC WellSpa iO device or ageLOC WellSpa iO Wireless Charger.
- Do not operate while bathing or showering. Do not submerge your device in water.
- To reduce the risk of shock, irritation, fire, or injury, do not place or store your device or charger where it may fall or be pulled into water, a bathtub, sink, or toilet.
- Do not modify the device in any way.

CHARGING YOUR AGELOC® WELLSPA iO™ DEVICE

The battery indicator will turn amber when the battery is low.

- Before charging, make sure the device and charger are completely dry.
- Plug the charger into a USB adapter and then into a wall outlet. Place the device on the charger to begin charging.

- While the device charges, the battery indicator breathes green.
- When the device is fully charged, the battery indicator will turn solid green.
- Do not keep the charger plugged in when not in use or after device has reached full charge.

NOTE

When the battery is critically low, the battery indicator will flash amber, and the device will turn off when the user attempts to turn it on. This indicates there is not enough battery remaining to perform a session.

If there are no indicator lights, charge the device. Indicator lights will turn on once the device is charged.

Our rechargeable battery helps reduce waste. We invite you to learn more about our sustainability initiatives by visiting nuskin.com/sustainability.

CONNECTING AGELOC® WELLSPA iO™ TO THE NU SKIN VERA® APP

To unlock the full ageLOC WellSpa iO experience, scan the QR code to download the Nu Skin Vera® app from the App Store or Google Play. Follow the app instructions to connect to your device.

PAIRING YOUR AGELOC® WELLSPA iO™

1. Download the Nu Skin Vera® app.
2. Ensure your phone's Bluetooth®* is enabled.
3. Sign in to Nu Skin Vera.
 - Already have a Nu Skin account? Use it to sign in.
 - New to Nu Skin? Sign up—it's easy!
4. Tap "Devices" in the navigation bar.
5. Tap "Pair New Device"
6. If the device is off, press the Power button to turn the device on.
7. Tap "Connect"
8. The device will pair and you're ready to go.

* The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by NSE Products, Inc. is under license. Other trademarks and trade names are those of their respective owners.

App Store and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. Google Play is a trademark of Google LLC. For information about how Nu Skin collects, uses, and shares Personal Data please visit <https://www.nuskin.com/content/global-privacy.html>

USING YOUR AGELOC® WELLSPA iO™ DEVICE SYSTEM

You may use your ageLOC WellSpa iO device system up to 5 minutes, three times per week to daily, on each body area.

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Smooth your serum or gel liberally on the targeted area, avoiding face, neck, chest, and groin. If you plan to target more than one area, only apply the product to the area to be immediately targeted.

STEP 2. STARTING A SESSION

Press the power button to turn on the device. All five indicator lights will illuminate the first time you power on. The number of illuminated indicator lights indicates the session time. Each time you press the Power button, you select a new time. Press the Power button again until you have selected your preferred time from 1–5 minutes. Your device will remember your setting for the next session. Place the device on the area of your skin where you've applied the serum or gel. The device will emit a tone, and your session will start automatically. Begin a new session for each targeted area.

STEP 3. DURING THE SESSION

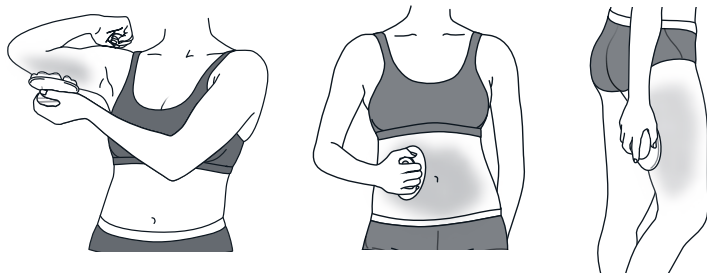
Lightly glide the device over your skin. Take care to avoid sensitive areas. Each minute will be counted down by the sound of a tone and an indicator light turning off.

- A lack of product or a loss of contact with your skin will interrupt the current causing WellSpa iO to pause and emit a tone. Resume by placing the metallic nodes on your skin and/or applying more product as needed to ensure a smooth glide.
- The device will emit an alert tone if you are moving the device too quickly.
- To turn the device off at any time, simply press and hold the Power button until the lights turn off.

STEP 4. SESSION COMPLETION

At the end of your session, the device will emit a tone and automatically turn off. You will hear the same sound when you manually power off your device. Massage any excess serum or gel into your skin.

See diagrams below for suggested usage areas and Nu Skin Vera® for full usage instructions.

CLEANING AND CARING FOR YOUR AGELOC® WELLSPA iO™

Use a slightly damp, lint-free wipe(s) to wipe away any serum, gel, or debris from the handheld device and charger after each use. Continue to wipe until visibly clean. Store in a cool, dry place. Do not immerse device in water or any liquid at any time. Using dry, non-linting wipe(s), thoroughly dry the external surfaces of the handheld device and charger.

TROUBLESHOOTING

LOSS OF CONDUCTIVITY DURING USE

If the metallic nodes lose contact with your skin, the current will be interrupted, causing WellSpa iO™ to pause, emit a tone, and the indicator lights to breathe.

PRO TIP: A lack of serum or gel may cause a loss of conductivity. Apply more product to your skin and ensure the device is in contact with both your skin and the product for optimal use.

MOVEMENT ALERTS DURING USE

The device will emit a sound if you are moving too quickly. The session will automatically resume when you resume the correct movement speed. Slowly glide the device over the targeted areas.

RESET BLUETOOTH® CONNECTION

If you have problems pairing with Nu Skin Vera® or if you need to unpair your device from your phone, conduct a Bluetooth reset.

1. Disconnect from Nu Skin Vera within the app by navigating to Device Settings > Forget this Device.

2. Disconnect your phone.

- **iPhone:** Navigate to your phone's Settings > Bluetooth > My Devices > NuSkinBDiO > Forget This Device.
- **Android:** Navigate to your phone's Settings > Connections > Bluetooth > Settings icon > Forget/Unpair.

3. Place WellSpa iO on the charger.

4. Press and hold the Power button for 5 seconds.

Wait for:

- Indicator lights 1-5 to flash once.

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by NSE Products Inc. is under license. Other trademarks and trade names are those of their respective owners.

FACTORY RESET

To factory reset your ageLOC® WellSpa iO™:

1. Disconnect from Nu Skin Vera® within the app by navigating to Device Settings > Forget this Device.
2. Disconnect from your phone.
 - **iPhone:** Navigate to your phone's Settings > Bluetooth® > My Devices > NuSkinBDiO > Forget This Device.
 - **Android:** Navigate to your phone's Settings > Connections > Bluetooth® > Settings icon > Forget/Unpair
3. Place WellSpa iO on the charger.
4. Press and hold the power button for 10 seconds.

Wait for:

- Indicator lights 1–5 to flash once.

Keep pressing and holding until:

- All indicator lights flash 3 times (at 10 seconds).

NOTE

ageLOC WellSpa iO automatically stores quality and usage information. When the device is factory reset, some device usage data will be retained for quality purposes. To view Nu Skin's privacy notice, visit:

-https://www.nuskin.com/en_US/corporate/privacy.html

DEVICE AND BATTERY DISPOSAL

You must dispose of ageLOC® WellSpa iO properly according to local laws and regulations. The device contains electronic components and a lithium-ion battery, which requires disposal separate from household waste. When ageLOC WellSpa iO reaches its end of life, contact local waste management authorities to learn about proper disposal and recycling options.

Dispose of batteries properly for your area.

MAINTENANCE

See the “Cleaning and Caring for Your ageLOC® WellSpa iO” section for routine maintenance.

Do not ship damaged devices by air.

SERVICING

Do not attempt to repair the device yourself, as this will void any warranty. There are no user-serviceable parts inside.

Please see the warranty section for service.

REPLACEMENT AND WARRANTY INFORMATION

Nu Skin guarantees your device to be free from defects in materials and workmanship for a period of two years from the original date of delivery to you. This warranty does not cover damage to the product as a result of misuse or accident. If the product becomes defective within the two-year warranty period, please call your local Nu Skin Support to arrange for a repair or replacement. You may be required to present a formal sales receipt showing the delivery date. Nu Skin reserves the right to inspect the device. This warranty does not affect your statutory rights in law specific to the country of purchase, such rights remain protected. This warranty is valid in Austria, Belgium, Czech Republic, Denmark, Finland, France, Germany, Hungary, Iceland, Ireland, Italy, Luxembourg, Norway, Poland, Portugal, Romania, Slovakia, South Africa, Spain, Sweden, Switzerland, The Netherlands, and United Kingdom.

PATENTS

US and international patents pending.

TECHNICAL AND REGULATORY INFORMATION

ageLOC® WellSpa iO™

Maximum Output:

22V, 640 μ A, 0.414Hz

 IPX5  IPX7

ageLOC WellSpa iO Wireless Charger

Input: 5V  min 500mA

 IPX4

DEFINITION OF SYMBOLS

Symbol	Definition
	Warning
	Dispose of Batteries Properly for Your Area
	Two Levels of Shock Protection
	Charge only with specified charger
	Direct Current

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© 2022 NED Product Europe BV. Alle Rechte vorbehalten. NED Product Europe BV ist ein eingetragenes Unternehmen in Belgien. NED Product Europe BV ist ein eingetragenes Unternehmen in Belgien. NED Product Europe BV ist ein eingetragenes Unternehmen in Belgien.
ZA Nu Skin Enterprises South Africa (PTY) LTD, 22 Magwa Crescent,
UK Nu Skin United Kingdom LTD, 11th Floor, 3 Piccadilly Place, Manchester,
©2023 NSE Products, Inc • www.nuskin.com