



NU SKIN® GALVANIC SPA® SYSTEM II

BRING SPA RESULTS HOME

NOW FEATURING  ageLOC™

BRING SPA RESULTS **home**

With the **stress of modern life**, spa treatments have become more than just occasional pampering, they are a **requisite** for **looking and feeling** great. The Nu Skin® Galvanic Spa® System II brings the day spa home to you. Using **cutting-edge**, patent-pending technology with the **Galvanic Spa**, **no appointment** is necessary to **smooth** the appearance of lines and wrinkles, **rejuvenate** your complexion, **revitalize** your scalp, and **renovate** your body.

For more than 50 years, spa professionals have used galvanic current in customized treatments to refresh and energize the skin. Through a gentle massaging action, this current helps focus cellular energy and enhance circulation.

Now you can get these same beneficial treatments in the convenience of your own home with the Nu Skin Galvanic Spa System II.

THE GALVANIC DIFFERENCE

- This patent-pending, hand-held device uses low level galvanic (electric) current to facilitate the transport of key ingredients into the skin for improved cleansing, hydration, circulation, and radiance.
- Galvanic Spa® Treatment Gel now features new ageLOC technology, a proprietary ingredient blend that can help slow a newly discovered source of accelerated aging in the epidermis—for results you can see now with benefits that last into the future.
- The gentle massaging action of each of the four attachments (face, treatment, body, scalp), combined with galvanic current, enhances circulation and draws vital nutrients to the outer layers of the skin.
- Studies have shown that a Nu Skin® Galvanic Spa® II treatment may enhance active ingredient delivery for up to 24 hours.





TREATED
(Results after two weeks
of treatment, three times a week.)

UNTREATED



TREATED
(Results after five weeks
of treatment, three times a week.)

UNTREATED



TREATED
(Results after four weeks
of treatment, three times a week.)

UNTREATED

REJUVENATE YOUR COMPLEXION

WITH GALVANIC SPA® FACIAL GELS WITH ageLOC

These **rejuvenating facial gels** work with the face conductor to draw out dulling debris and allow **pores to breathe**, while delivering **beneficial nutrients** that focus cellular energy and facilitate recovery from stress.

The revolutionary **ageLOC** ingredient blend in Treatment Gel also helps **slow the visible signs of aging** at their **source**.

Your skin will look **visibly younger**

and feel **vibrant** and **revitalized**.

ageLOC™



- **PRE-TREAT FACIAL GEL:**
 - A cleansing agent purifies the skin and allows pores to breathe.
 - Marigold, sea kelp, and vanilla extracts soothe the skin.
- **TREATMENT FACIAL GEL:**
 - Proprietary ageLOC blend helps slow free radical production at its source in the epidermis and inhibits the visible signs of aging.*
 - Magnesium promotes cellular energy.
 - Arginine facilitates recovery from stress.
 - Aspartate naturally balances your skin's pH.
 - Tocopheryl Nicotinate, a form of vitamin E, helps moisturize, condition, and provide antioxidant protection.

* See page 6 for further information.

LET OUR SCIENCE BE YOUR SECRET



Time will never tell.

average appearance of **7** years older

High arNOX →

+18
+16
+14
+12
+10
+8
+7
+6
+4
+2
Actual age
-2
-4
-6
-7
-8
-10
-12
-14
-16

Low arNOX →

average appearance of **7** years younger


ONLY NU SKIN'S EXCLUSIVE, PATENT-PENDING AGELOC TECHNOLOGY CAN HELP TO SLOW FREE RADICAL PRODUCTION AT ITS SOURCE IN THE EPIDERMIS—HELPING YOU CONTROL YOUR SKIN'S FUTURE.

Present on the surface of our cells, an invisible enemy—arNOX—is capable of generating skin damaging free radicals around the clock that can also destroy collagen and elastin and diminish production of these critical proteins. The result? Looking years older than you are.

A recent study by Stanford University has shown that individuals with low arNOX levels are clinically graded to appear to look an average of seven years younger than their actual chronological age. Individuals with high levels of arNOX appear to be an average of seven years older than their actual age.


THE SCIENCE BEHIND NU SKIN® GALVANIC SPA® II

GALVANIC ACTION PRE-TREAT

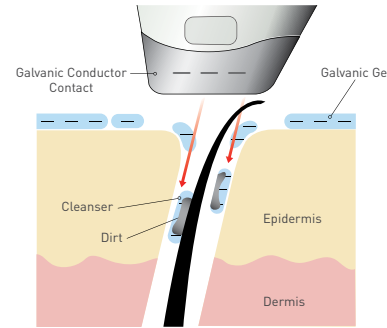

- During the **Pre-Treat Cycle**, the Galvanic Instrument and the Pre-Treat Gel are **negatively charged**.
- These negative charges repel each other, **delivering the Pre-Treat Gel into the skin**.
- The Pre-Treat Gel binds to impurities in the skin.

GALVANIC ACTION TREATMENT

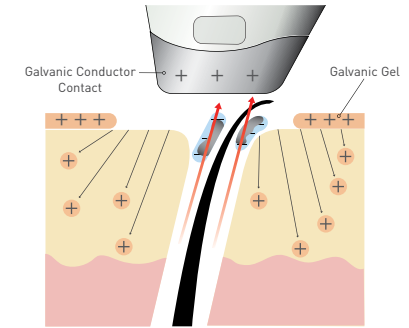

- During the **Treatment Cycle**, the Galvanic Instrument and the Treatment Gel are **positively charged**.
- These positive charges repel each other helping to **deliver the beneficial ingredients to the skin**, including the exclusive ageLOC ingredient blend.
- The positively charged instrument attracts the remaining negatively charged impurities from the Pre-Treat cycle, drawing them out of the skin.

FACIAL TECHNIQUES


ITS' ALL IN THE MOTION

With all of these movements, gently come back to the starting position and repeat several times.

- ① EYE AREA: Following the bone structure around your eye, start under the eye and move up around the eye, holding firmly just under the eyebrow.
- ② SMOOTH SMILE LINE: Start just under the nose, move to the smile line and hold and then move straight upward and outward to the temple and hold.
- ③ FIRM JOWLS: Start in the middle of the chin, below your lips, move up to the jowl line and hold and then continue in an upward and outward motion to the edge of your face just in front of your ear.
- ④ SMOOTH FOREHEAD: Start at the bridge of your nose, move the instrument up to your “thinker lines” and hold, then continue on to your hair line. Move across the forehead using these same movements.
- ⑤ CALM PUFFINESS: Hold the instrument at a 90-degree angle on the side of your nose and then move to the middle of your cheekbone and back.
- ⑥ SMOOTH NECK: Start at the base of your throat. Move upward to the jaw line, hold for a few seconds, glide the unit back down to the starting position and then move upward and outward, holding each time at the jaw line, and then returning to the hollow of the throat.

REMEMBER

- Keep your hand moist, with your fingers always touching the chrome panel.
- When you touch the instrument to your skin you will hear one to three beeps as the instrument sets the current to your skin's needs.
- Always keep the instrument in contact with your skin. If you lose contact, return the instrument to your skin and resume.
- Keep the facial conductor flat against the skin.
- Work the instrument in upward, outward movements.

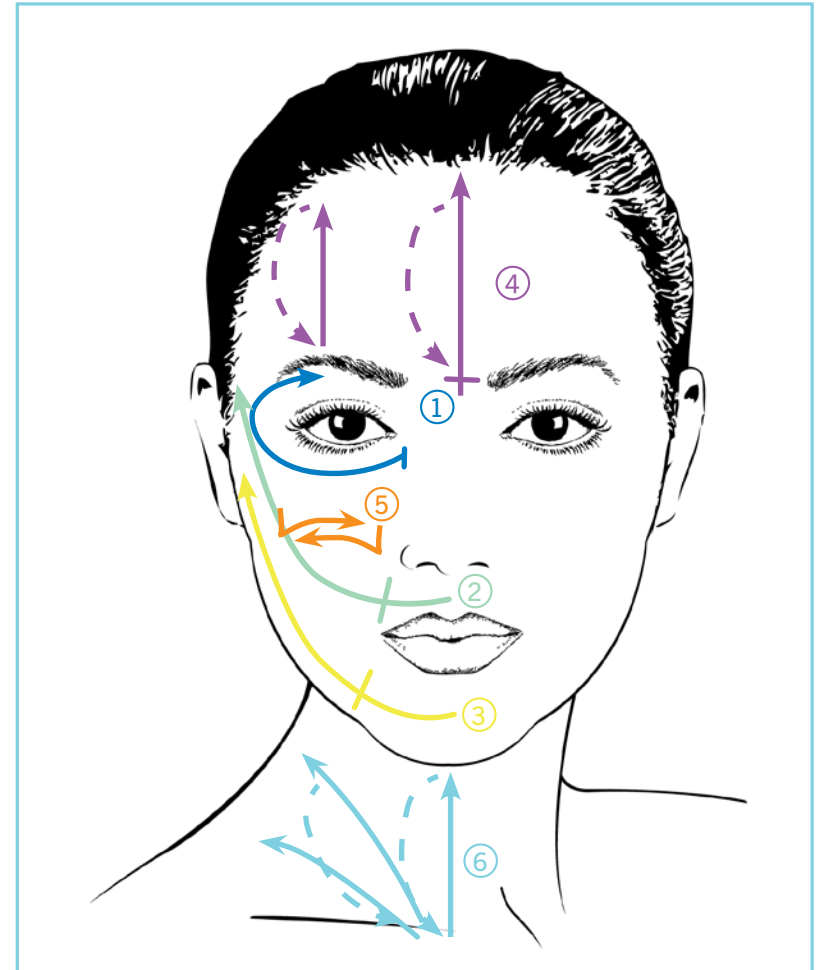

FOCUS ON YOUR AREAS OF CONCERN.

NU SKIN® GALVANIC SPA® II + ACTIVE SKIN CARE INGREDIENTS =

MORE EFFECTIVE ABSORPTION


Research using the Nu Skin Galvanic Spa System II provides exciting new information. A recent study* by a third-party lab measured the ability of the Galvanic Spa to help skin care ingredients delivery to the skin.

Once the galvanic treatment was concluded, measurements were taken at three different intervals. The first measurement was taken at **0 hours**, the second at **6 hours**, and the last was taken at **24 hours** after completing a galvanic treatment.

Active "A" is an active skin care ingredient found in Nu Skin treatment products. You can see that at both **6 and 24 hours** there is an **increase** in the ability of Active "A" to enter the galvanic-treated skin. Galvanic treatment increases absorption making active ingredients more available for skin improvement.

To capitalize on these benefits, the Galvanic Spa® should be used two to three times a week.

SKIN CARE INGREDIENT DELIVERY WITH GALVANIC TREATMENT

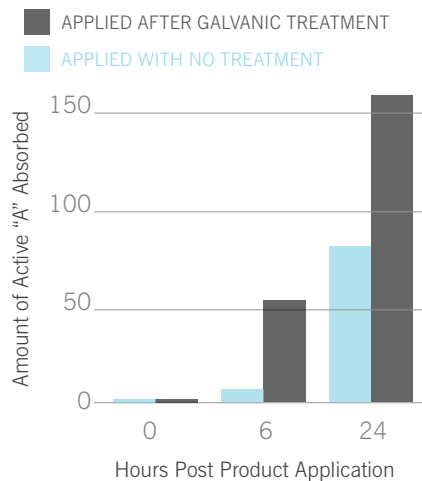

*The results show the impact of using quality skin care in combination with the Galvanic Spa® System II. Coupled with the Nu Skin 180® Anti-Aging Skin Therapy System, the Galvanic Spa can dramatically improve the look and feel of your skin.

Research studies conducted at the University of Oklahoma, 2006


ARE YOU MAXIMIZING YOUR GALVANIC SPA® II BENEFITS?

GET YOUR BOUNCE BACK WITH TRU FACE® ESSENCE ULTRA



Tru Face Essence Ultra is a daily use product that helps to **BUILD** and **PROTECT** skin's elastin content. The key ingredient, Ethocyn®, is a unique anti-aging molecule that helps improve the skin's elastin content, which is critical in maintaining smooth, firm, youthful looking skin. Ethocyn has been clinically proven in a third-party study to improve elastin fiber content to levels typically found in 18–25 year olds.

ELASTIN CONTENT OF OLDER SKIN TREATED WITH ETHOCYN VS. UNTREATED YOUNG SKIN ($p < 0.001$)

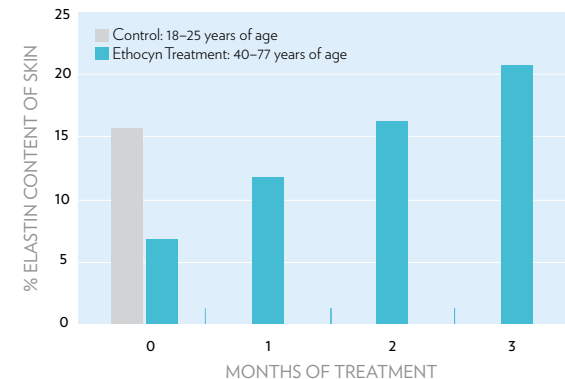


SMOOTH APPEARANCE OF LINES AND WRINKLES

WITH TRU FACE® LINE CORRECTOR

Galvanic Spa® II with the new Spot Treatment conductor works synergistically with **Tru Face Line Corrector**, helping to **enhance the delivery** of a pro-collagen peptide. Peptides help **reduce the appearance** of moderate to deep **lines** by sending age-reversing signals to collagen-producing cells. The **new spot treatment conductor** allows you to better **target specific areas** of the face and gives an **added boost** to your existing **anti-aging** treatment regimen.

- PALMITOYL PENTAPEPTIDE-3 (PRO-COLLAGEN PEPTIDE)
 - Used with the Galvanic Spa II instrument, this peptide is more effectively absorbed into the skin, helping to deliver enhanced results.
 - Sends age-reversing signals to collagen-producing cells to help reduce the appearance of moderate to deep lines.
- JUST AS A KEY CAUSES THE ENGINE OF A CAR TO START, PEPTIDES ARE THE KEY TO TURNING ON COLLAGEN PRODUCTION.
- COLLAGEN STRUCTURE STRENGTHENS THE SKIN, REVERSING SIGNS OF AGING AND PROMOTING YOUNGER LOOKING SKIN.

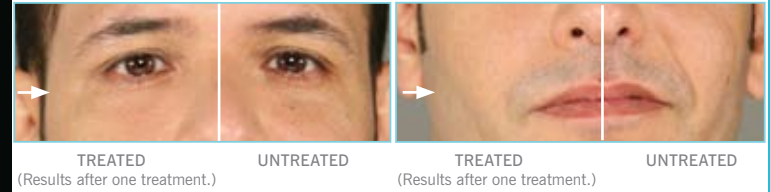

In a recent study, dermatological grading showed that using the Galvanic Spa Instrument with Tru Face Line Corrector improved the appearance of the following key areas:

- Fine Lines: 20%
- Crow's Feet: 21%
- Tautness: 17%



UNTREATED TREATED
(Results after one treatment.)



UNTREATED TREATED
(Results after one treatment.)

RENOVATE YOUR BODY

WITH GALVANIC SPA® II BODY SHAPING GEL



This **exclusive body gel** works together with the ergonomically designed body conductor to **target pockets of fat** and **toxins**, minimizing and **smoothing the appearance of cellulite** for a firmer, more **toned body**.

- THEOBROMINE (THEOBROMA COCOA EXTRACT)
– Contains a molecule shown to target enzymes found in fat cells. Targeting these enzymes may help activate fat breakdown.
- GOLDEN CHAMOMILE (CHRYSANTHELLUM INDICUM)
– Rich in flavonoids shown in testing to inhibit the storage of fat.
- RESEARCH SHOWS THAT WHEN GALVANIC SPA II BODY SHAPING GEL IS USED TWICE A DAY, THE BENEFITS ARE CUMULATIVE.




UNTREATED ARM

TREATED ARM
(Results after four weeks of treatment, three times a week.)

- Smooth Appearance of Cellulite
- Firm Skin




Use DERMATIC EFFECTS® Body Contouring Lotion daily to extend results between your Galvanic Spa® II Body Shaping Gel treatments.

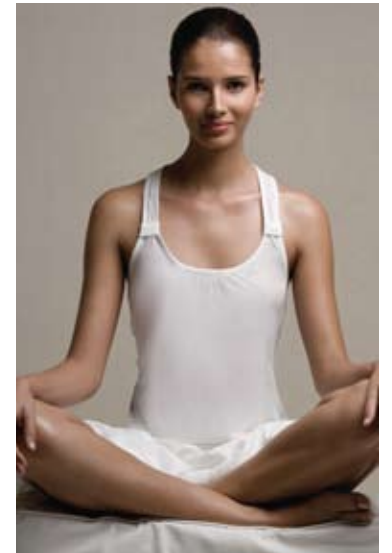

REVITALIZE YOUR SCALP AND HAIR

WITH NUTRIOL® HAIR FITNESS TREATMENT

This **technologically advanced** treatment helps your hair look and feel **more abundant**. Nutriol Hair Fitness Treatment helps facilitate the transport of key ingredients into the scalp, promoting **healthy** and **abundant hair**.

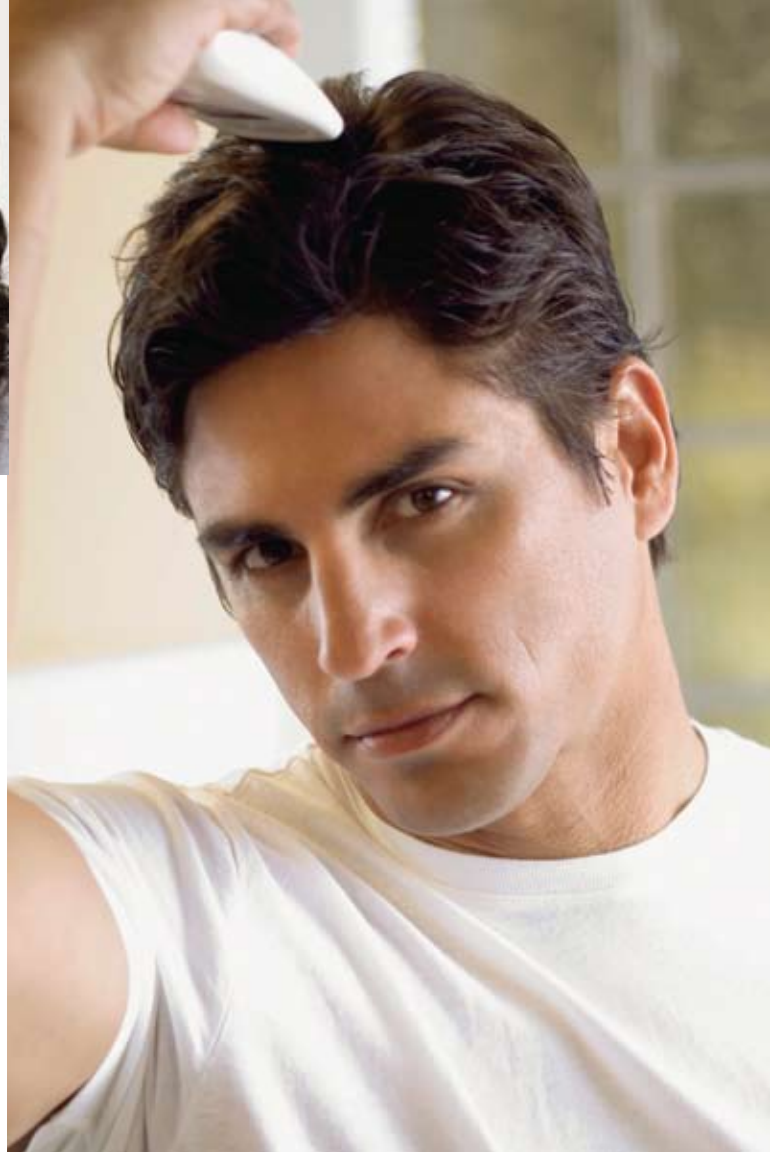

- TRICALGOXYL®
 - Rich in kelp sulphated oligosaccharides and is a clinically proven patented scalp-nourishing ingredient.
 - Works with other essential ingredients to hydrate the hair and remineralize the scalp, helping hair regain its natural health, strength, and shine—for maximum vitality and abundance.



Nutriol® Shampoo contains the active ingredient Tricalgoxyl® and works to help hair look and feel more abundant between galvanic treatments.

Specially formulated products work synergistically with each proprietary attachment to deliver targeted benefits. This combination provides remarkable health and beauty benefits for the face, scalp, and body.



Galvanic Spa® System II

Professional spa-like treatments—no appointment necessary. Rejuvenate your complexion, revitalize your scalp, and renovate your body with the patented Galvanic Spa System II.

WHITE US 01003867 \$300
CA 02003867 \$330

BLACK US 01101319 \$300
CA 02101319 \$330



Galvanic Spa® Facial Gels with ageLOC™

Slow visible signs of aging at the source. Newly reformulated Facial Gels contain all of the great benefits of the original facial gels but now Treatment Gel also features ageLOC technology, a proprietary new ingredient blend that helps slow the signs of aging at their source. Get results that you can see now with benefits that last into the future.

US 01003876 \$46 4 vials each
CA 02100403 \$50.60 4 vials each



Galvanic Spa® II Body Shaping Gel

Take action and smooth away the appearance of cellulite. Featuring theobromine and golden chamomile, this advanced gel helps promote a firmer, more toned appearance. For optimal results, use with the Galvanic Spa II Instrument.

US 01102725 \$44 5 oz.
CA 02102725 \$48.75 150 mL



Nutriol® Shampoo

Head-turning results for your head of hair. Featuring Tricalgoxyl®, this daily shampoo is formulated to restore fullness and shine.

US 01101219 \$33 4.2 oz.
CA 02101219 \$36.15 125 mL



Nutriol® Hair Fitness Treatment

Helps hair look and feel more abundant. This leave-in treatment remineralizes the scalp and hydrates the hair with Tricalgoxyl®, a clinically proven seaweed derivative.

US 01101218 \$63 12 vials
CA 02101218 \$69.20 12 vials



Dermatic Effects® Body Contouring Lotion

The perfect partner for fighting cellulite. An eight-week clinical study showed a 55% reduction in the appearance of cellulite, a 37% improvement in skin firmness, and a 69% enhancement in skin smoothness.

US 01110270 \$44 5 oz.
CA 02110270 \$48.75 150 mL



Tru Face® Line Corrector

Smooths lines with the power of peptides. This advanced formula features a revolutionary pro-collagen peptide that helps minimize the appearance of moderate to deep lines around your mouth, eyes, and forehead in just four weeks.

US 01102704 \$48 1 oz.
CA 02102704 \$53.50 30 mL

ORDER NOW

1.800.487.1000 US
1.800.877.6100 CA

nuskinusa.com
nuskin.com/canada

ageLOC™ SPA PACKAGE

A special collection of Galvanic products. See nuskinusa.com for details.

US Item 01 103433 White
US Item 01 103477 Black

CA Item 02 103433 White
CA Item 02 103477 Black



nuskinusa.com 1.800.487.1000
nuskin.com/canada 1.800.877.6100

5 pack 01 006833/4 USA
5 pack 02 006833/4 CAN

 **NU SKIN**[®]

©2008 Nu Skin International, Inc. • Provo, UT 84601